ENCHEN

Mint 3

Portable Oral Irrigator

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
ENCHEN

Mint 3

Portable oral irrigator · three intensity levels · USB-C charging

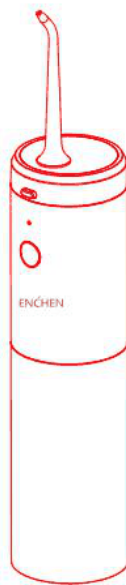
User Manual

Read this manual carefully before using the product.

Contents

1	Product Description	3
2	Use	4
3	Important Warnings	6
4	Care and Maintenance	7
5	Faults and Their Remedies	8
6	Technical Data	9
7	Hazardous Substance Content	10

1. Product Description

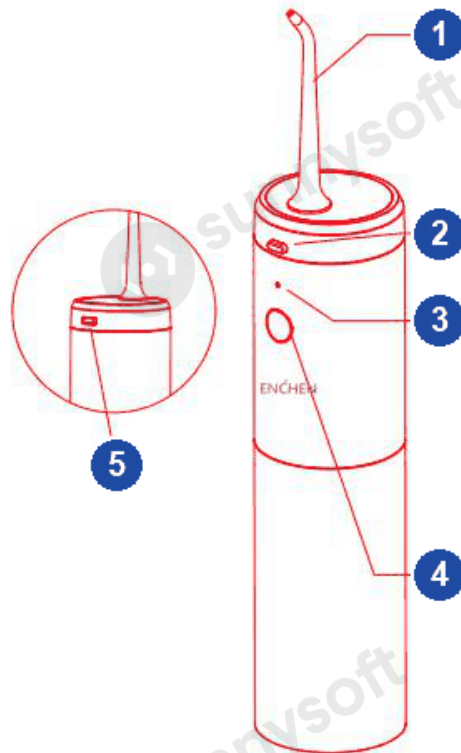
- 1 nozzle
- 2 nozzle release button
- 3 power indicator
- 4 power button / level selector
- 5 USB-C charging port (on the back of the product, in the detail on the left)

What the Indicator Shows


During normal operation the indicator lights up white. If the battery is charged to less than 30 %, it lights up yellow. While charging, the white indicator flashes; once the product is fully charged, it stays white for a longer period.

2. Use

1. Push the nozzle into its opening until you hear a click. Turn the water tank anticlockwise to remove it and fill it. Tap water is recommended; if it is too cold, use warm water below 40 °C. Do NOT add salt, toothpaste or any other chemicals to the tank. Fit the filled tank back on by turning it clockwise. When fitting it, do NOT overtighten the tank or turn it upside down, so that the water does not spill.

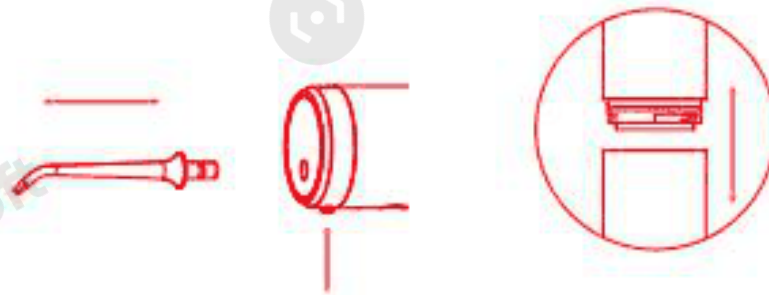2. While using the product, lean slightly forward, open your mouth a little and let the water flow out freely, or spit it out. While rinsing, do NOT stand upright or sit down, so that the water does not splash onto your face and body. For the best result, start at the molars and work towards the front teeth, tooth by tooth along the gum line. Pay attention to the pits and fissures of the teeth, to the wisdom teeth and to other hard-to-reach places. As you move the nozzle along the gum line, pause briefly at every interdental space so that the area between the teeth is cleaned as well.

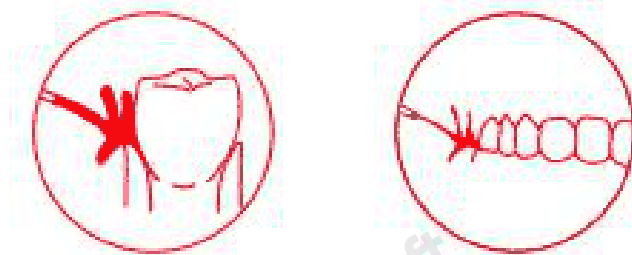3. Press the power button briefly to switch between the first, second and third level and switching off.

First level (gentle)

After you press the button, the product starts at this level and the operating indicator lights up white.

Second level (standard)

Press the button again within 3 seconds of starting the first level. The product switches to this level and the operating indicator lights up white.

Third level (strong)


Press the button again within 3 seconds of starting the second level. The product switches to this level and the operating indicator lights up white.

Note


If you press the button after more than 3 seconds of operation at any level, the product switches off. To preserve battery life, the product switches itself off by a protective function after 3 minutes of continuous operation.

4. After use, let the product run for 2 to 3 seconds so that the remaining water drains out. Open the cover of the charging port, insert one end of the charging cable into the port of the product and the other end into the power adapter, then plug the adapter into a socket.

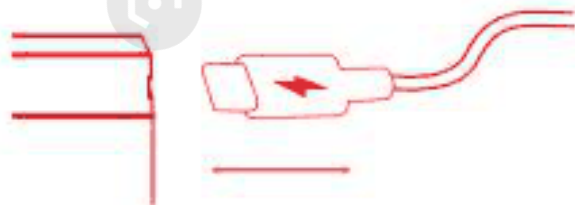Internal Layout

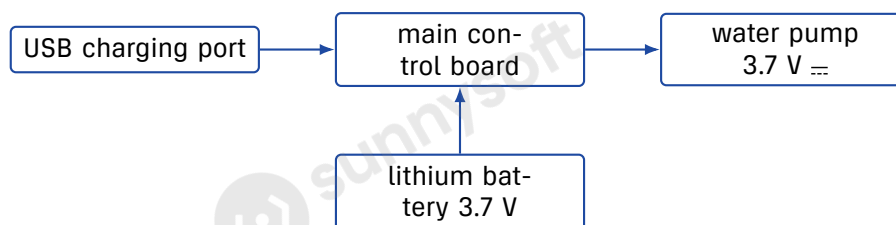3. Important Warnings

To reduce the risk of failure, electric shock, injury, fire, death or damage to the product and property, observe the following safety precautions.

1. Do not disassemble the product or modify its parts at will – you could damage it or cause an electric shock. Unauthorised disassembly voids the warranty.
2. Do not put water hotter than 40 °C into the tank, to avoid scalding.
3. If the product falls into water while charging, first unplug the power adapter from the socket and only then take the product out, to avoid an electric shock.
4. Do NOT throw the product into fire and do not heat it. Do NOT charge it in a place with high temperature or humidity, do NOT leave it in such a place for a longer time and do not expose it to direct sunlight. While charging, do NOT pull out the power and charging plug with wet hands.
5. A power adapter is not included. For charging, use an adapter with a DC output of 5 V / 1 A that complies with the applicable standards. Before using the adapter, check that the stated voltage matches the local mains voltage.
6. If you do not use the product for a longer time, unplug the power and charging plug to prevent fire or other hazards.
7. Keep the product clean, use running or warm water to rinse your mouth, and after use let the product run for 2 to 3 seconds so that the water drains out.
8. The product has no replaceable or repairable parts. If a repair is necessary, return the product to the manufacturer.
9. Do not use the product without the nozzle fitted.
10. The product is suitable only for persons over 6 years of age. Children should be accompanied by an adult.
11. Charge the product at room temperature.
12. When charging at a low temperature (below 15 °C), the operating time may be shortened. If the temperature is too low (below 0 °C), charging may not be possible.
13. Use running tap water or warm water, and press the button only after you have brought the nozzle close to your teeth.
14. The built-in rechargeable battery may contain pollutants. When disposing of the product, remove the battery and hand it in at a designated collection point, or hand in the whole product at a designated collection point.

15. If the battery cannot be removed, you can send the product to the manufacturer's service centre, where they will help you remove it and dispose of it in an environmentally friendly way.
16. Do NOT remove the rechargeable battery at will. You may remove it only when disposing of the product.

4. Care and Maintenance

Warning

Do NOT immerse the product in water in order to clean it.

- Use clean water or a neutral detergent for cleaning. Do NOT use corrosive or abrasive agents, they could damage the product. Do NOT wash the product with hot water above 50 °C.
- Remove the nozzle after every use and clean the product.
- The whole device can be rinsed under running water and then stored away.
- Wipe the outer surface with a soft damp cloth to remove dirt, and after every use wipe the relevant parts dry with a clean dry cloth.
- Dentists recommend replacing the nozzle once a month.

5. Faults and Their Remedies

Fault	Possible cause	Remedy
The product does not work.	The product has just been purchased.	Connect the product to the power supply and charge it for 3 hours until it is fully charged.
	The product has not been used for several months.	Contact customer service.
After charging it runs for only a few minutes.	The battery is low.	Connect the product to the power supply and charge it for 3 hours until it is fully charged.
	The battery has reached the end of its life.	Contact customer service.
Low water pressure.	Dirt is stuck in the nozzle.	Replace the nozzle.
	There is dirt in the water pump.	Contact customer service.
No water comes out.	The tank has run out of water.	Fill the tank with water.
	There is dirt in the water pump.	Contact customer service.

If you have further questions, contact customer service.

6. Technical Data

Product name	portable oral irrigator
Model	Mint 3
Rated voltage	3.7 V $\equiv$
Rated power input	4 W
Power adapter	not included; 5 V / 1 A $\equiv$

The conformity marks in the image are taken from the manufacturer's documentation.

Manufacturer

Shanghai Qiaopusi Intelligent Technology Development Co., Ltd.

Room 103, Building 3, No. 320 Donglan Rd, Minhang District, Shanghai, China, 201102

7. Hazardous Substance Content

The table was compiled by the manufacturer in accordance with regulation SJ/T 11364.

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○ means that the content of the substance in all homogeneous materials of the part is below the limit set by standard GB/T 26572. × means that the content of the substance in at least one of the homogeneous materials of the part is above this limit. The manufacturer states that the product meets the requirements of the RoHS Directive and that no advanced technology is yet available that would replace or reduce the lead content in electronic ceramics, optical glass, steel and copper alloys. The number on the product label indicates that the usual period of use from the environmental protection point of view is 10 years.

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